

C-FILL Endodontic Obturation system MINI


Precise temperature control
Perfect filling



Order name: C-FILL
Mini



The tip of plunger with minimum diameter 0.45mm, 02 taper, can press the effective filling position that is 3 mm from the Apex.

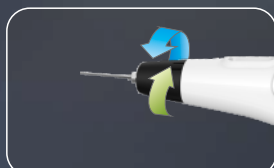

The gun head can be rotated by $\pm 100^\circ$ for easy use.



Double-sided button operation, convenient and fast



360° dynamic charging display, more humane.

Model: C-Fill

Order code: C-FILL Mini G type

Battery: DC3.7V 2200mAH

Heating time: 15 s To 200 °C.

Working temperature: 150 °C, 180 °C, 200 °C, 230 °C

Adapter: input AC 100~240V, output DC5V 1.5A

Use time 1.5H

Needle specifications: 23G, 25G



Accurate temperature control

Automatically cool down to 150°C on the base
Save electricity and prevent glue overflow

Endodontic Obturation Systems



Faster heating
G type can archive 200°C in 15s



New wrench design
Help you to push more accurately



Order code: C-FILL Mini P type

Battery: DC3.7V 2000mAH

Heating time: 4S(With heat protection function, the single heating time is up to 4 S)

Working temperature: 150 °C, 180 °C, 200 °C, 230 °C

Adapter: input AC 100~240V, output DC5V 1.5A

Use time 1.5H

Plunger specifications: F 45/02, FM 50/04, M 55/06, ML 55/08



Different angle
Easy to reach any treatment places



Faster heating
P type can archive 200°C in 0.5s

Intelligent reminder
Sound and flashing prompt
With protection function, single heating time up to 4s



Precise temperature control
Good performance give you the most smooth cutting of gutta-percha



C-FILL mini P type

Heat plunger		
Size	Color code	Taper/Size
#XF	Green	0.02#35
#F	White	0.02#45
#FM	Yellow	0.04#50
#M	Red	0.06#55
#ML	Red	0.08#55



Needle		
Gauge	Diameter	Length
23ga	0.65mm	23mm
23ga	0.65mm	23mm
25ga	0.55mm	23mm
25ga	0.55mm	23mm

