

THE DENTAL
SOLUTIONS
COMPANY™



Intego & Intego Pro

Performance you can rely on

dentsplysirona.com/intego





Your reliable partner

Intego is the best treatment center of its class. It is based on the innovative power and award-winning design of our existing treatment centers. As a reliable partner for an efficient practice, Intego has all the features that provide optimal support: ergonomic design for better workflows, high quality standards in processing and design, and a great value – this is why Intego is our performance class.

Great value

Excellent performance and innovative functions at a great price, customized options for great value.

Ergonomic design

Designed for intuitive handling and optimized, integrated workflows.

Quality made in Germany

Defined quality assurance processes, durable materials, and sturdy, future-proof design.

It's up to you

One size fits all – we know that in a dental practice, very often this is not true. This is why Intego comes in two versions to match your individual requirements.

Intego

If you are looking for a treatment center to optimally support your personal workflows at a fair price, Intego is your choice. While being simple to operate, it brings all relevant functions where you need them to be. The

ergonomic design offers high comfort to both you and your patients. The Compact water unit with a fresh water bottle makes it easy to comply with all infection prevention regulations.

Intego Pro

The highlight of Intego Pro undoubtedly is its comfort water unit, that brings a comprehensive infection prevention concept to your practice. The integrated disinfection system reduces the workload of your team

through automatic purging and sanitization of the water lines, as well as a permanent disinfection with Dentosept S.



Intego TS (OTP arm design)

- Ergonomic positioning of the dentist element due to a wide range support arm system
- Functional tray system with detachable tray that can be positioned flexibly



Great value you can rely on

Your work day is complex and places many different demands on the treatment center you choose. Intego gives you the perfect solution: Choose from two treatment concepts that you can configure flexibly. Whether you select a whip arm or a hanging hoses design - every option offers optimal performance at a fair price.

CS - The whip arm design



- Ideal for the widest range of treatment situations with four- or two-handed dentistry
- Can be rotated up to 240° for optimal positioning
- Maximum distance between the dentist element and instruments is 90cm

- 1 Sprayvit E**
Optimal visibility due to integrated LED light.
- 2 BL ISO E motor**
Sterilizable brushless motor with speed and torque control during treatment. Especially lightweight, compact and with a wide speed range.
- 3 SiroSonic L scaler**
Handpiece and tips can be sterilized. Controlled, straight oscillation with electronic power and water control.
- 4 SiroCam F or AF or AF+**
Detailed intraoral images with the easy to operate SiroCam F or AF/AF+ with autofocus.



Intuitive operation

Difficult treatments require your full concentration. This is why Intego is intuitive – with two state-of-the-art user interfaces and innovative instruments for stress-free treatment. So you can concentrate on what's really important – your patient and a perfect treatment result.

State-of-the-art, intuitive interfaces



EasyPad

Clear and intuitive user interface – standard in Intego. Basic Endo function available with EasyPad.




EasyTouch

This user interface can be optionally integrated. More operating comfort thanks to intuitive touchscreen technology and additional operating and display options. Required for advanced Endo function.

Integration possibilities for endodontic treatments

The integrated endodontic function of Intego is the ideal solution for dental professionals looking for more safety, efficiency, and simplicity in performing endodontic procedures. Having everything integrated results in a streamlined and absolutely hygienic workflow and reduces set up times to a minimum.

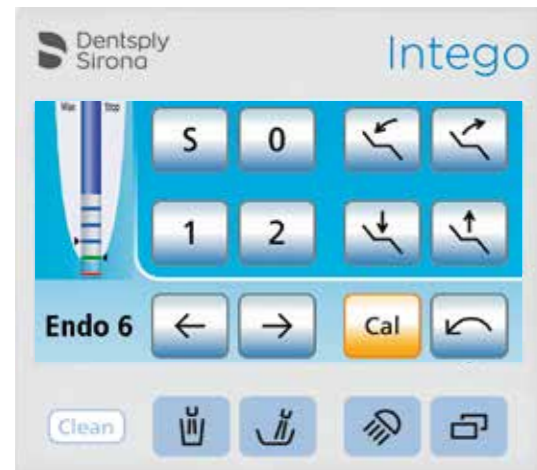


ApexLocator

An ApexLocator can also be integrated into the EasyTouch to allow for permanent control of the working depth during endodontic treatment and a reduction of X-ray exposure.




Heliodent Plus unit model

Fast, reliable intraoral X-rays with the Heliodent Plus unit model. The durable X-ray unit gives you maximum flexibility and optimal image quality.

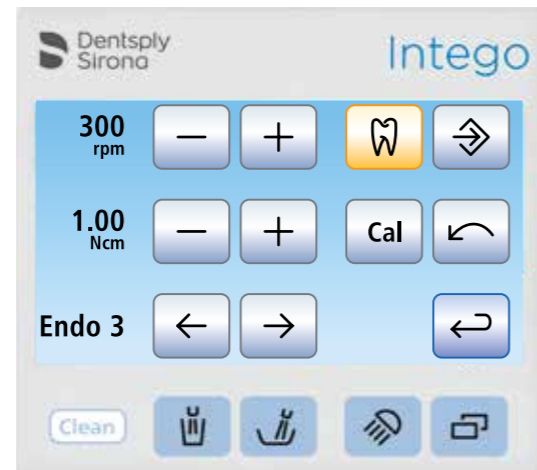

Integrated Torque Control

Speed and torque can be preset for 6 different files, and the motor stops automatically when the set maximum torque is reached. This ensures greater safety for endodontic treatments.

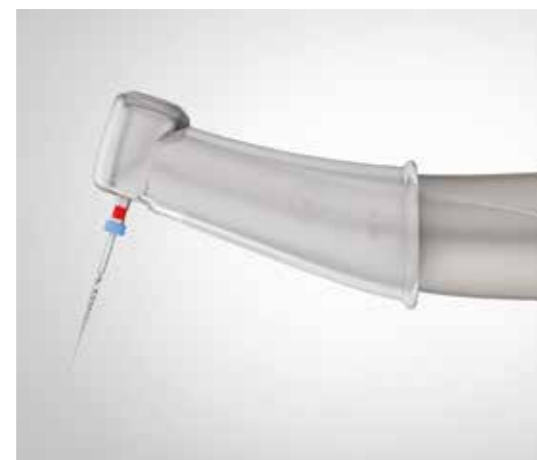

T1 Line contra-angle handpiece Endo 6 L with light

Improved visibility during endodontic treatment increases the rate of a successful treatment. Another available option without light is the contra-angle handpiece Endo 6:1.

Reciprocal file systems

Intego integrates the Dentsply Sirona Maillefer and VDW reciprocal file systems, enabling a fast and easy use of common file systems.

WaveOne®

WaveOne® Gold
Reciprocating Files

Reciproc®
One file endo



Impressive patient comfort

Your patients will also appreciate this comfort: The Intego patient chair has an ergonomic design and, together with the headrest, ensures that your patients are especially comfortable. The chair lowers far enough for even shorter patients so that they can get in and out easily. And there is no doubt – patients who are comfortable feel more secure and remain loyal to your practice.



Dentsply Sirona's unique ErgoMotion for Intego

Your benefits

- Better positioning for maxillar or mandibular treatment
- Less overstretching of patient's head
- No overstretching and compressing of patient's spine
- No readjustment of headrest necessary

Double articulating headrest

Quickly adjustable to various treatment situations as well as patients with varying statures.

Premium/Lounge upholstery

The high-quality upholstery ensures comfort evenduring long treatments.

ErgoMotion

The tilting of the backrest matches the movement of the patient's spine to prevent stretching and compression effects.

Armrests

Your patients can rest their arms comfortably.

Lifting system

The unique chair lift drive gently moves patients into the proper position.

Vario footrest

The Vario footrest that comes with the Premium upholstery can be easily extended to provide even very tall patients with ergonomic positioning.

Safety you can rely on – with Intego

Dentists all over the world have the highest demands for the safety of their patients and practice team. The Intego infection prevention concept is designed so you can comply quickly and easily with all infection control regulations. The assistant element and water unit designs make infection prevention procedures even easier and less time consuming. All of this guarantees perfect infection control.

Intego with
compact water unit

Cuspidor

Removable and swiveling, easy to clean, and easily accessible.

Purge function

Manual purging of all instrument waterlines on the dentist and assistant side.

Suction systems

Compatible with all conventional wet and dry suction systems, including amalgam separator.

Intego

The Compact water unit and assistant element can be adapted to many different demands and requirements by freely selectable suction systems and infection prevention options.



Fresh water bottle

For controlled water quality. Can also be filled with Dentosept S solution for permanent disinfection and manual sanitization.



Automatic suction hose cleaning

Cleaning adapter for automatic flushing of saliva ejector and HVE.



Smooth surfaces

The largely gap-free design of the surfaces ensures maximum infection control in the shortest time possible.

Safety you can rely on – with Intego Pro

Dentists all over the world share the highest demands for the safety of their patients and practice team. The Intego Pro infection prevention concept is designed so you can comply quickly and easily with all infection control regulations. The assistant element and water unit designs make infection prevention procedures even easier and less time consuming. All of this guarantees perfect infection control.

Intego Pro with
comfort water unit

Cuspidor
Removable and swiveling, made
of high-quality glass.

**Purge and
Autopurge function**
The instrument hoses can
be flushed for a short time
or cleaned intensively via
the cleaning adapters.

Suction systems
Compatible with all conventional
central wet or dry suction systems,
including amalgam separator.

Intego Pro

The Comfort water unit and assistant element are based on a comprehensive infection prevention concept.



Integrated disinfection system
Permanent disinfection of all water
lines with Dentosept S. Automatic
sanitization function for regular
intensive disinfection. Can also be
used as a self-contained water
supply.



Integrated cleaning adapters
For automatic purging of
instrument hoses (Autopurge) and
sanitization of the water lines.



Chemical suction hose cleaning
Cleaning agent container for
chemical flushing of the suction
lines. The concentration of cleaning
agents can be individually adjusted.

Intego Ambidextrous – a masterpiece of flexibility

Ambidextrous option available for both
Intego and Intego Pro

Perfect ergonomic turning

In organizing an efficient everyday practice routine, it is essential that the practice equipment can be adjusted not only to your requirements, but also to the requirements of your team. This is an even greater challenge in a group practice, where a treatment center will be used by several dentists. The treatment area must be quickly and easily adjustable. This process must also be absolutely safe and stable. The Intego Ambidextrous is designed precisely for this purpose.

The turning process in a nutshell



1 The turn process is activated by a single step on the respective foot pedal and the backrest moves automatically in the designated Turn position.




2 As the water unit is smoothly moved behind the backrest on a robust support arm, the dentist element and OP light can be simultaneously swiveled to the new position with no additional space required.



3 Finally, if they are existent, tray holder and camera support can be swapped.

Your benefits

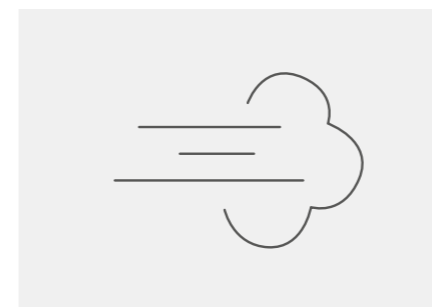


Fast

Due to the intelligent mechanism design, the treatment position can be switched in just 15 seconds.



Easy

There is no need for special training, the process can be intuitively handled by your team. The swivel mechanism is released with just one click on the foot pedal and the backrest travels automatically into the correct position.

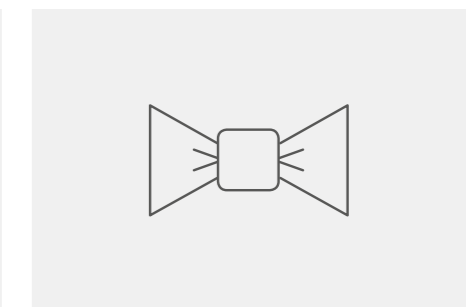


Elegant


The treatment center can be converted smoothly and effortlessly in a single flow. No jerky movements, no screws, no assembly.

More freedom of movement with no extra space needed

Intego Ambidextrous provides great flexibility in a small area. Despite the conversion option, no more space is required than for our existing Dentsply Sirona treatment centers. The intelligently designed rotating mechanism enables compact and space-saving movement at all times in every phase. During the conversion process, the treatment center remains completely stable in all positions. The water unit is guided safely and in a controlled manner around the backrest - this rules out the risk of damage. The following figures exemplify the fact that Intego Ambidextrous in every phase of the conversion movement does not require more space than other Intego treatment centers.

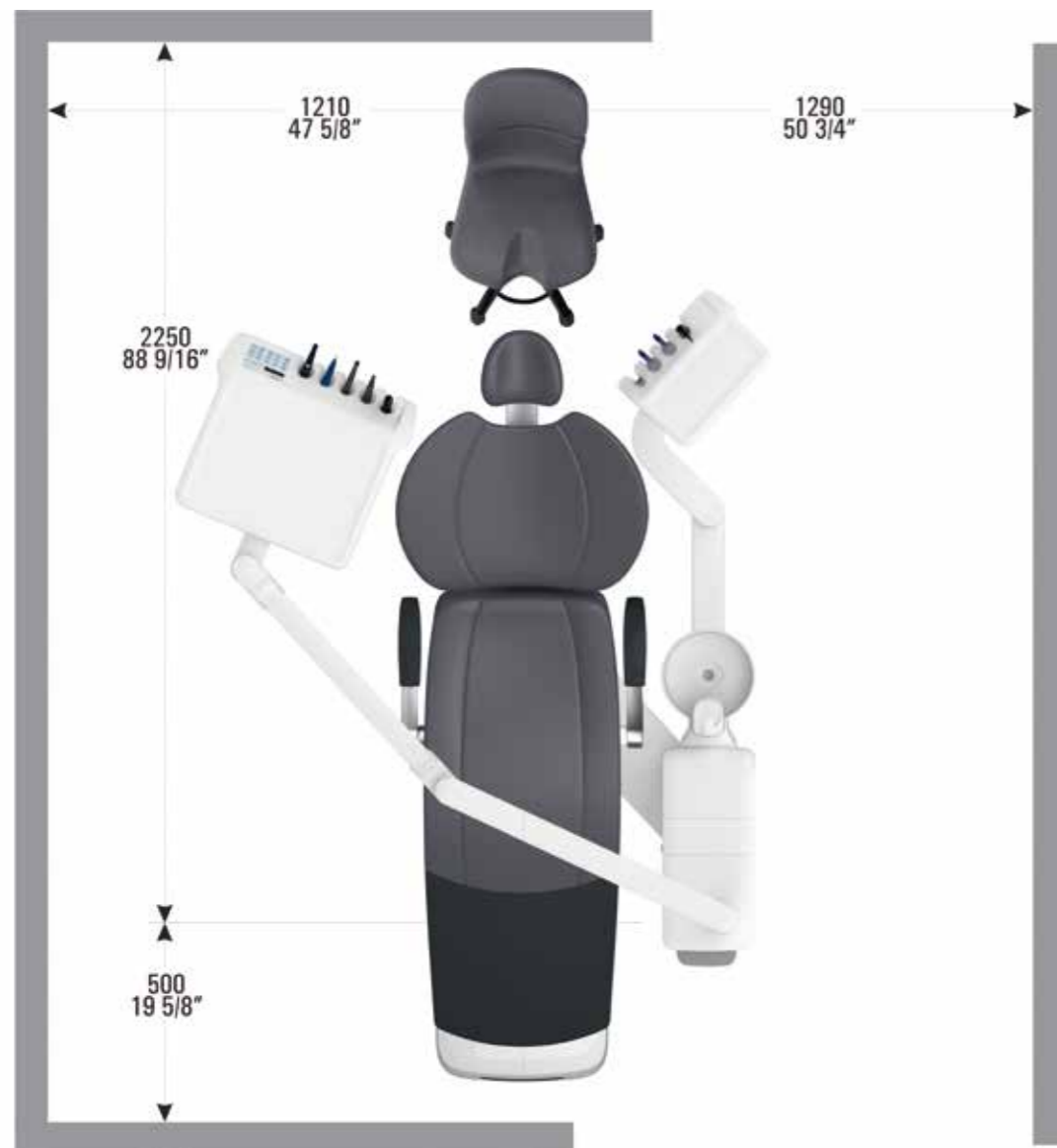

This figure illustrates the required space of any Intego when used in 12 o'clock treatment position.

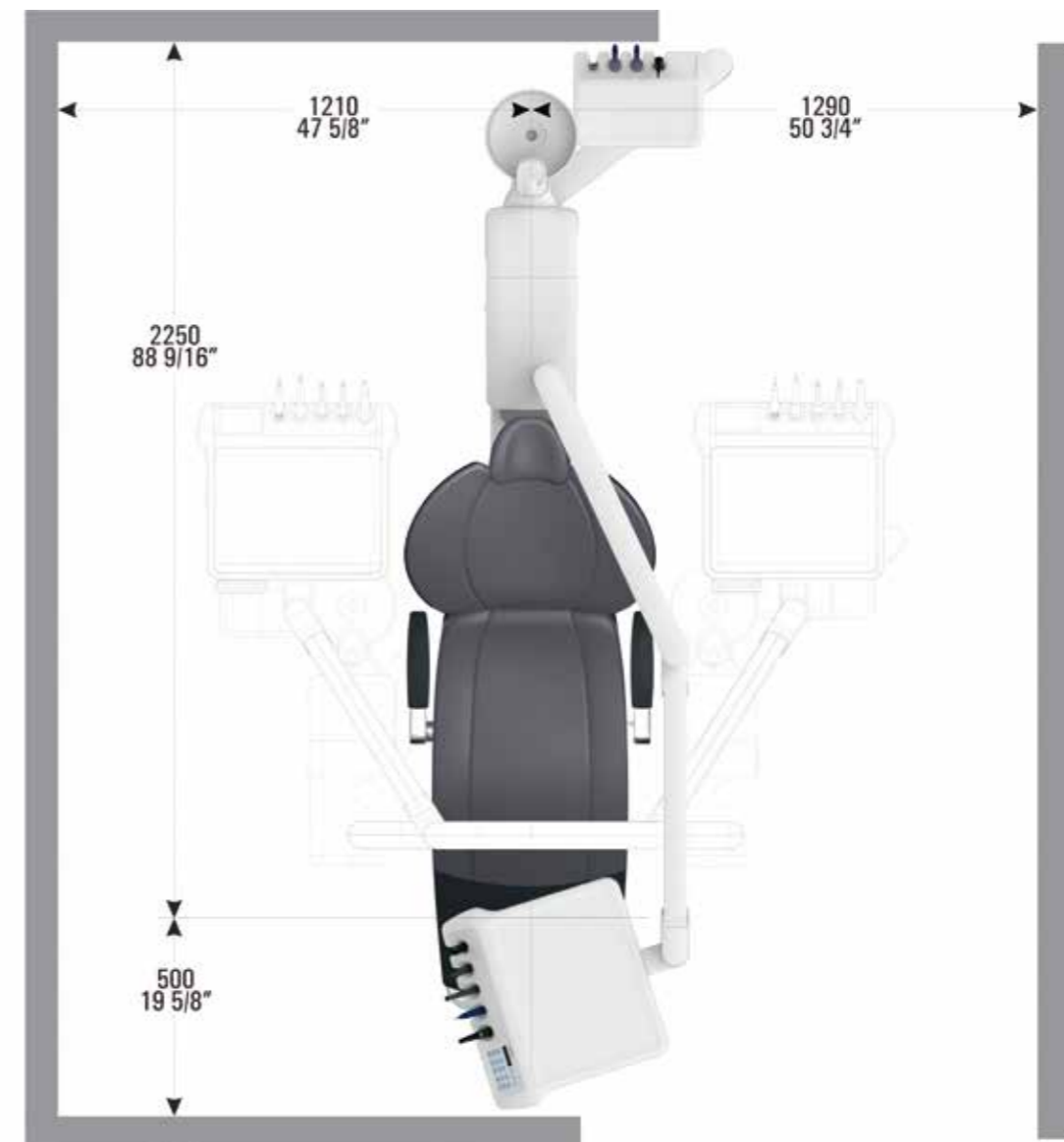

This figure illustrates the maximum required space of Intego Ambidextrous with Comfort assistant element during the conversion process.

Perfect ergonomic position for every treatment

Achieve perfect results even at the end of a long working day – with your expertise and the support of our treatment center. No matter if you prefer to treat your patients as a right- or left-handed practitioner, in every situation you will benefit from intuitive sitting, comfortable positioning, optimal visibility and integrated workflows. This combination will lead to perfect clinical results.



More freedom of movement

- Thin backrest for ideal freedom of leg movement
- Flexible positioning of the foot control
- Easy height adjustment via dentist element or optional foot switch
 - Swiveling cuspidor

Optimized instrument placement

- Perfect coordination of dentist and assistant element
- Flexible positioning of the additional tray for more storage place
- Compact and comfort assistant element based on your treatment concept

Comfortable working stools

- Hugo, Carl, Paul
- New: Saddle stool Theo for improved lower jaw treatment
- Supports healthy sitting position throughout the day

Optimum visibility

- Double-articulating headrest for quick and flexible adaptations
- Flat headrest with variable positions through the magnetic cushion
- LEDview Plus or LEDlight Plus operating light, LED instrument light

Design meets function

Prestigious design awards such as the iF communication award and the red dot design award bear witness to the outstanding design quality of the Dentsply Sirona Treatment Centers. Intego offers an enormous range of design options for your practice with three different color schemes. Let your imagination run free and put together the perfect combination.



Upholstery colors

The elegant world



Classic, timeless combinations for a truly sophisticated look.

The natural world



Soft, warm shades for a soothing treatment atmosphere.

The vital world



Vivid colors and bright tones for a modern, lively atmosphere.

* Available colors for Lounge upholstery.

Stitching colors



Standard

Red

Dark grey

Light blue

Rose

Upholstery options



Premium upholstery

The high-quality, seamless upholstery ensures comfort even during long treatments.



Lounge upholstery

Soft, pleasant surface material and comfortable filling for the utmost relaxation. Contoured surface for perfect positioning.

Flexible into the future

In a good practice, customized settings add to your comfort. This is why you can configure the Intego Pro exactly the way you need for optimal support of your treatment workflow. But flexibility means even more to us – and you can prepare Intego Pro with technical options and interfaces so you are ready for future applications.

Many models

Choice of two headrests

The flat headrest can be adjusted quickly and flexibly. Optionally, you can select a double articulating headrest that is easier to position and allows you optimal visibility even in difficult-to-see areas.



Flat headrest



Double articulating headrest

Two operating light options – always optimal illumination

Two choices of LED lights for optimal illumination of the operating field. Our premium operating light LEDview Plus offers with its multi color LEDs and sophisticated technical features natural light in perfection.



LEDlight Plus



LEDview Plus

Assistant element: Comfort or Compact

Both assistant elements are very space-saving and ensure maximum freedom of movement. With the Comfort model, you can work comfortably even without an assistant.



Comfort assistant element



Compact assistant element

Ultramodern equipment

The SiroCam F/AF/AF+ autofocus cameras combine superior image quality with easy operation.



	SiroCam F	SiroCam AF	SiroCam AF+
Focus	Fixed Focus	Auto Focus	Auto Focus
Image format 16:10	✓	✓	✓
Pixels (usable)	0.42 megapixel	0.65 megapixel	1.00 megapixel
Workspace	✓ 5–20 mm	✓ 2–130 mm	✓ 2–130 mm
Resolution live-image	824 x 514	824 x 514	824 x 514
Resolution still-image	824 x 514	1024 x 640	1276 x 796

Future-proof design

Interfaces: flexible into the future

A network interface in the patient chair and a USB port in the dentist or assistant element allow you to integrate new functions and network your practice – troubleshooting by your service technician is also much easier.



Integrated USB interface

Sivision 22" AC monitor

With this HD monitor, you can explain intra-oral camera and X-ray images, software and planning views, videos and PowerPoint presentations to your patient directly at the treatment center.



Sivision 22" AC monitor

Quality you can rely on

We don't want to just satisfy our customers, we want to impress them. This is why quality, performance, and service are key for us. Dentsply Sirona treatment units combine high quality with real innovations that perceptibly enhance treatment comfort and workflow.

Reliability and durability



Precise, elaborate quality tests and robust materials make Intego Pro a reliable treatment center with an above-average service life.

High quality



Intego Pro immediately shows our attention to detail. The materials used – robust, shapeable plastics, solid steel and aluminum components on the inside, high-quality glass and upholstery – meet the high quality demands of a premium product.

High standards



Dentsply Sirona is a technology leader with more than 130 years of innovation history. Our engineers and scientists in Research & Development not only pay close attention to innovation, but also quality.

Future-proof design



At the Dentsply Sirona innovation center, future trends are anticipated and products are prepared to be equipped for the technology demands of the future.






Individual configuration

Intego

Basic unit
<p>Patient chair</p> <ul style="list-style-type: none"> Premium upholstery Pneumatic foot control Flat headrest Ambidextrous option <p>Dentist element</p> <ul style="list-style-type: none"> 5 instrument positions EasyPad user interface Supply for 2 drill drives <p>Assistant element</p> <ul style="list-style-type: none"> Compact version on a swivel arm directly mounted to the cuspidor 4 handpiece positions Saliva ejector and HVE <p>Water unit</p> <ul style="list-style-type: none"> Standard wet suction Ready for integration of a fresh water bottle Cuspidor with manual swivel function

Version
<p>TS version with hanging hoses</p> 
<p>CS version with whip arms</p> 

User packages
<p> Start</p> <ul style="list-style-type: none"> Standard 3-way syringe Light source in the dentist element 2 turbine hoses LEDlight Plus*
<p> Advance</p> <ul style="list-style-type: none"> C+ electronic foot control Sprayvit E in the dentist element Light source in the dentist element 1 turbine hose 1 BL ISO E micro-motor with supply electronics LEDlight Plus Upcharge option: EasyTouch
<p> Innovation</p> <ul style="list-style-type: none"> C+ electronic foot control Sprayvit E in the dentist element Light source in the dentist element 1 turbine hose 1 BL ISO E micro-motor with supply electronics SiroSonic L scaler LEDview Plus Upcharge option: EasyTouch

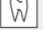Option packages
<p> Sivision Ready</p> <ul style="list-style-type: none"> USB interface for the dentist or assistant element Network interface
<p> Sivision Start with/without monitor</p> <ul style="list-style-type: none"> USB interface for the dentist or assistant element Network interface SiroCam F (SiroCam AF or AF+ at additional charge) Monitor mount on the light support column Optional: 22" AC monitor
<p> Ergo</p> <ul style="list-style-type: none"> ErgoMotion Chair lift system 185 kg 2 armrests Double articulated headrest
<p> Endo Pro**</p> <ul style="list-style-type: none"> Integrated torque control ApexLocator Endo 6:1 contra-angle handpiece e-Control
<p>* LEDview Plus at additional charge ** Only with EasyTouch</p>




Availability and options of basic unit and packages can vary by country.

Intego Pro

Basic unit
<p>Patient chair</p> <ul style="list-style-type: none"> Premium upholstery Pneumatic instrument foot control Flat headrest 4-way foot switch Ambidextrous option <p>Dentist element</p> <ul style="list-style-type: none"> 5 instrument positions EasyPad user interface Supply for 2 drill drives Instrument light source Supply electronics for up to 2 motors and scaler <p>Assistant element</p> <ul style="list-style-type: none"> Comfort assistant element 4 handpiece positions Saliva ejector and HVE Light source for Sprayvit E <p>Water unit</p> <ul style="list-style-type: none"> Ready for the integration of an automatic disinfection device Integrated cleaning adapters Cuspidor with manual swivel function Wet suction with multi-station valve Integrated water flush suction hose cleaning Water heater

Version
<p>TS version with hanging hoses</p> 
<p>CS version with whip arms</p> 

User packages
<p> Start Pro</p> <ul style="list-style-type: none"> C+ electronic foot control Sprayvit E in the dentist element 1 turbine hose 1 BL ISO E micro-motor LEDlight Plus
<p> Advance Pro</p> <ul style="list-style-type: none"> C+ electronic foot control Sprayvit E in the dentist element 1 turbine hose 1 BL ISO E micro-motor SiroSonic L scaler EasyTouch user interface LEDview Plus
<p> Innovation Pro</p> <ul style="list-style-type: none"> C+ electronic foot control Sprayvit E in the dentist element Third drill drive 1 turbine hose 2 BL ISO E micro-motors SiroSonic L scaler EasyTouch user interface Sprayvit E in the assistant element LEDview Plus

Option packages
<p> Sivision Ready</p> <ul style="list-style-type: none"> USB interface for the dentist or assistant element Network interface
<p> Sivision Start with/without monitor</p> <ul style="list-style-type: none"> USB interface for the dentist or assistant element Network interface SiroCam F (SiroCam AF or AF+ at additional charge) Monitor mount on the light support column Optional: 22" AC monitor
<p> Endo Pro**</p> <ul style="list-style-type: none"> Integrated torque control ApexLocator Endo 6:1 contra-angle handpiece e-Control
<p>* LEDview Plus at additional charge ** Only with EasyTouch</p>

Availability and options of basic unit and packages can vary by country.

Your individual options

Here you can view the complete list of individual options, and a checklist of what is included with basic units.

Intego

Patient chair		Ambi-dextrous
2 treatment positions, S and O position	■	■
Integrated junction box	■	■
Chair lift system 140 kg	■	■
Chair lift system 185 kg	□Er	□Er
Flat headrest	■	■
Double articulating headrest	□Er	□Er
Premium upholstery	■	■
Lounge upholstery	□	□
ErgoMotion	□Er	□Er
4-way foot switch	□	□
Booster cushion for small patients	□	□
Right arm rest	□Er	□Er
Left arm rest	□Er	□Er
Valve set for external device	□	□
Vacuum pressure reduction kit	□	□
Cover installation kit	□	□
Separate junction box	□	□
Adapter plate	□	□

Dentist element		Ambi-dextrous
Height-adjustable support arm with pneumatic brake	■	■
5 instrument positions	■	■
Additional camera holder for TS dentist element	□	–
Hydro-pneumatic supply for 2 drill drives w/o fiber optics	■	■
Pneumatic foot control	■	■
C+ electronic foot control	□A, I	□A, I
EasyPad user interface with 2 user-specific settings	■	■
EasyTouch user interface with 4 user-specific settings	□T	□T
Sprayvit E with light	□A, I	□A, I
Sprayvit E with light (water button on the left)	□A, I	□A, I
Standard 3-way syringe	□S	□S
Turbine hose with light	□Sx2, A, I	□Sx2, A, I
Borden turbine hose	□	□
1 st BL ISO E motor with hose	□A, I	□A, I
2 nd BL ISO E motor with hose	□	□
SiroSonic L scaler	□I	□I
Third drill drive	□	□
Intra-oral light source kit for turbine and Sprayvit E	□S, A, I	□S, A, I
Torque control for endodontic treatments	□E	□E
eControl: Speed range of 100-40.000 rpm	□E	□E
Integrated ApexLocator	□E	□E
Endo 6:1 contra-angle handpiece	□E	□E
Intra-oral x-ray viewer	□	□
Purge and manual sanitization function	□	□
Flexible TS tray holder for 2 standard size trays	–	□
Additional TS tray holder for 2 standard size trays	□	–
Silicone pad for TS dentist element	□	□
Additional CS tray holder for 1 standard size tray	□	□
Additional CS tray holder for 2 standard size trays	□	□
Silicon pad for CS dentist element	□	□
Removable silicone sleeves for handles	□	□

Assistant element		Ambi-dextrous
Compact assistant element, swivel mounted directly on the cuspidor arm	■	■
Comfort assistant element, ideal for two-handed dentistry	□	□
4 instrument positions	■	■
Saliva ejector and HVE	■	■
Sprayvit E with light	□	□
Sprayvit E with light (water button on the left)	□	□
Standard 3-way syringe	□	□
2 nd HVE kit	□	□
Satelec Mini LED curing light	□	□
Touch pad	□	□
Light source for Sprayvit E	□	□
Silicone pad for Comfort assistant element	□	□

Water unit		Ambi-dextrous
Prepared for fresh water bottle	■	■
Cuspidor with manual swivel function	■	■
White cuspidor	■	■
Water unit without cuspidor	□	□
Prepared for standard wet suction	■	■
Prepared for standard wet suction w/ automatic suction control	□	□
Multi-station valve	□	□
Cuspidor valve for standard wet suction	□	□
Automatic separator	□	□
Amalgam separator	□	□
Fresh water bottle w/ol switch to public water supply	□	□
Fresh water bottle with switch to public water supply	□	□
Water heater	□	□
Water flush suction hose cleaning system	□	□

Operating light/Heliodont Plus options		Ambi-dextrous
LEDlight Plus	□S, A	□S, A
LEDview Plus	□I	□I
LEDview Plus ceiling light	□	–
Adaptation for LEDview Plus/Heliodont Plus ceiling combination	□	–
Adaptation for Heliodont Plus unit model	□	–

Patient communication		Ambi-dextrous
USB interface in the dentist or assistant element	□SiR, SiS	□SiR, SiS
SiroCam F	□SiS	□SiS
SiroCam AF	□	□
SiroCam AF+	□	□
Monitor mount on the light support column	□SiS	□SiS
22" AC monitor (with integrated power supply unit)	□SiS*	□SiS*
Network interface	□SiR, SiS	□SiR, SiS

■ Included in basic unit
 □ Option
 – Not possible/applicable
 SiR Sivation Ready package
 SiS Sivation Start package
 SiS* Sivation Start package with monitor

S Start Pro package
 A Advance Pro package
 I Innovation Pro package
 E Endo Pro package
 Er Ergo package

Availability of options and features can vary by country.

Intego Pro

Patient chair		Pro Ambi-dextrous
2 treatment positions, S and O position	■	■
Integrated junction box	■	■
Chair lift system 185 kg	■	■
Flat headrest	■	■
Double articulating headrest	□	□
Premium upholstery	■	■
Lounge upholstery	□	□
ErgoMotion	□	□
4-way foot switch	■	■
Booster cushion for small patients	□	□
Right arm rest	□	□
Left arm rest	□	□
Valve set for external device	□	□
Vacuum pressure reduction kit	□	□
Cover installation kit	□	□
Separate junction box	□	□
Adapter plate	□	□

Dentist element		Pro Ambi-dextrous
Height-adjustable support arm with pneumatic brake	■	■
5 instrument positions	■	■
Additional camera holder for TS dentist element	□	–
Flexible TS camera holder	–	□
Hydro-pneumatic supply for 2 drill drives w/o fiber optics	■	■
Pneumatic foot control	■	■
C+ electronic foot control	□S, A, I	□S, A, I
EasyPad user interface with 2 user-specific settings	■	■
EasyTouch user interface with 4 user-specific settings	□A, I	□A, I
Sprayvit E with light	□S, A, I	□S, A, I
Sprayvit E with light (water button on the left)	□S, A, I	□S, A, I
Turbine hose with light	□S, A, I	□S, A, I
1 st BL ISO E motor with hose	□S, A, I	□S, A, I
2 nd BL ISO E motor with hose	□I	□I
SiroSonic L scaler	□A, I	□A, I
Third drill drive	□I	□I
Satelec Mini LED curing light	□	□
Intra-oral light source kit for turbine and Sprayvit E	■	■
Torque control for endodontic treatments	□E	□E
eControl: Speed range of 100-40.000 rpm	□E	□E
Integrated ApexLocator	□E	□E
Endo 6:1 contra-angle handpiece	□E	□E
Intra-oral x-ray viewer	□	□
Purge and Auto Purge	■	■
Flexible TS tray holder for 2 standard size trays	–	□
Additional TS tray holder for 2 standard size trays	□	–
Silicone pad for TS dentist element	□	□
Additional CS tray holder for 1 standard size tray	□	□
Additional CS tray holder for 2 standard size trays	□	□
Silicon pad for CS dentist element	□	□
Removable silicone sleeves for handles	□	□

Availability of options and features can vary by country.

Assistant element		Pro Ambi-dextrous
Compact assistant element, swivel mounted directly on the cuspidor arm	–	■
Comfort assistant element, ideal for two-handed dentistry	■	□
4 instrument positions	■	■
Saliva ejector and HVE	■	■
Sprayvit E with light	□I	□I
Sprayvit E with light (water button on the left)	□I	□I
2 nd HVE kit	□	□
Satelec Mini LED curing light	□	□
Touch pad	□	□
Light source for Sprayvit E	■	■
Silicone pad for Comfort assistant element	□	□

Water unit		Pro Ambi-dextrous
Prepared for integrated disinfection device	■	■
Cuspidor with manual swivel function	■	■
White cuspidor	■	■
Glass cuspidor	□	–
Water unit without cuspidor	□	□
Multi-station valve	■	■
Cuspidor valve for standard wet suction	□	□
Automatic separator	□	□
Amalgam separator	□	□
Integrated disinfection system with permanent water-line disinfection and automatic sanitization function	□	□
Water heater	■	■
Water flush suction hose cleaning system	■	■
Chemical suction hose cleaning system	□	□

Operating light/Heliodont Plus options		Pro Ambi-dextrous
LEDlight Plus	□S	□S
LEDview Plus	□I, T	□I, T
LEDview Plus ceiling light	□	–
Adaptation for LEDview Plus/Heliodont Plus ceiling combination	□	–
Adaptation for Heliodont Plus unit model	□	–

Patient communication		Pro Ambi-dextrous
USB interface in the dentist or assistant element	□SiR, SiS	□SiR, SiS
SiroCam F	□SiS	□SiS
SiroCam AF	□	□
SiroCam AF+	□	□
Monitor mount on the light support column	□SiS	□SiS
22" AC monitor (with integrated power supply unit)	□SiS*	□SiS*
Network interface	□SiR, SiS	□SiR, SiS

■ Included in basic unit
 □ Option
 – Not possible/applicable
 SiR Sivation Ready package
 SiS Sivation Start package
 SiS* Sivation Start package with monitor

S Start Pro package
 A Advance Pro package
 I Innovation Pro package
 E Endo Pro package

You decide what fits to your needs

This is an overview of numerous options that are available for Intego and Intego Pro. Your treatment center – your choice.



Hugo



Carl



Paul



C+ electronic foot control



LEDview Plus



LEDlight Plus



Theo



Additional tray holder CS



Additional tray holder TS



Valve set for external devices



Heliodent Plus unit model



22" AC monitor



Double articulating headrest



Flat headrest




Additional camera holder TS



SiroCam F



SiroCam AF

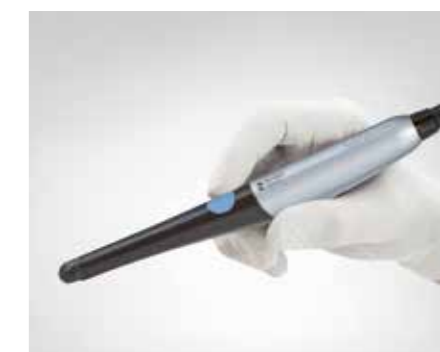

SiroCam AF+



X-ray viewer



Comfort assistant element



Compact assistant element



Saletec Mini LED curing light



Intego without cuspidor



4-way foot switch

Dentsply Sirona

Sirona Dental Systems GmbH
Fabrikstraße 31, 64625 Bensheim, Deutschland
dentsplysirona.com

Procedural Solutions

Preventive
Restorative
Orthodontics
Endodontics
Implants
Prosthetics

Enabling Technologies

CAD/CAM
Imaging
Treatment Centers
Instruments